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
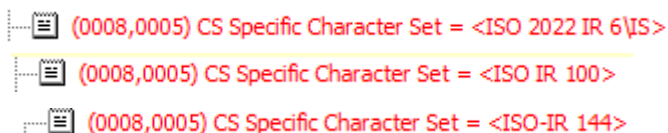
(0008,0005) CS Specific Character Set = <ISO 2022 IR 6\IS>
(0008,0005) CS Specific Character Set = <ISO IR 100>
(0008,0005) CS Specific Character Set = <ISO-IR 144>

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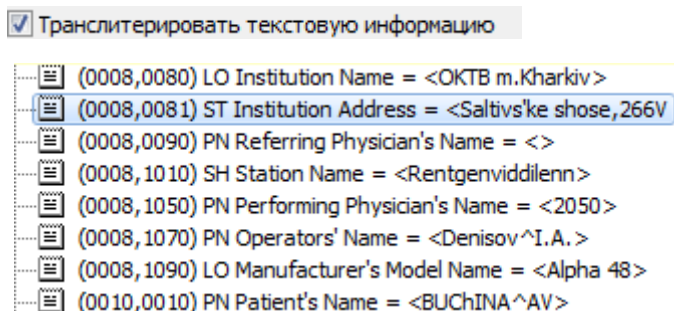
(0008,0080) LO Institution Name = <ОКТБ м.Харків>
(0008,0081) ST Institution Address = <Saltivs'ke shose,266V>
(0008,0090) PN Referring Physician's Name = <>
(0008,1010) SH Station Name = <Rentgenviddilenn>
(0008,1050) PN Performing Physician's Name = <2050>
(0008,1070) PN Operators' Name = <Denisov^I.A.>
(0008,1090) LO Manufacturer's Model Name = <Alpha 48>
(0010,0010) PN Patient's Name = <BUChINA^AV>

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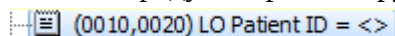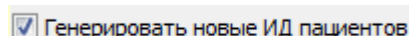
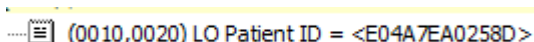
Новые Instance UID: Study Series SOP

 (0020,000D) UI Study Instance UID = <0>
 (0020,000E) UI Series Instance UID = <0>

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Error! Tags 0002,0003 and 0008,0018 are different, it's standard violation

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 (0020,000D) UI Study Instance UID = <1.2.276.0.7230010.3.1.2.21158767.2676.1473409496.603>
 (0020,000E) UI Series Instance UID = <1.2.276.0.7230010.3.1.3.21158767.2676.1473409496.604>

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COPYRIGHT

Unless otherwise specified, the DCMТK software package has the following copyright:

```
/*
 * Copyright (C) 1994-2011, OFFIS e.V.
 * All rights reserved.
 *
 * This software and supporting documentation were developed by
 *
 * OFFIS e.V.
 * R&D Division Health
 * Escherweg 2
 * 26121 Oldenburg, Germany
 *
 * Redistribution and use in source and binary forms, with or without
 * modification, are permitted provided that the following conditions
 * are met:
 *
 * - Redistributions of source code must retain the above copyright
 *   notice, this list of conditions and the following disclaimer.
 *
 * - Redistributions in binary form must reproduce the above copyright
```

```

* notice, this list of conditions and the following disclaimer in the
* documentation and/or other materials provided with the distribution.
*
* - Neither the name of OFFIS nor the names of its contributors may be
* used to endorse or promote products derived from this software
* without specific prior written permission.
*
* THIS SOFTWARE IS PROVIDED BY THE COPYRIGHT HOLDERS AND CONTRIBUTORS
* "AS IS" AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, BUT NOT
* LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR
* A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE COPYRIGHT
* HOLDER OR CONTRIBUTORS BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL,
* SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT
* LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE,
* DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY
* THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT
* (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE
* OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.
*
*/

```

Please note that some DCMTK modules, especially those that are not part of the free toolkit, are covered by a separate license which can be found in the COPYRIGHT file in the corresponding module directory.

Some portions of the DCMTK software package are derived from earlier versions of this software with the following copyright, and can be identified by the following copyright notice located in each source file:

```

/*
* Copyright (C) 1993/1994, OFFIS, Oldenburg University and CERIUUM
*
* This software and supporting documentation were
* developed by
*
* Institut OFFIS
* Bereich Kommunikationssysteme
* Westerstr. 10-12
* 26121 Oldenburg, Germany
*
* Fachbereich Informatik
* Abteilung Prozessinformatik
* Carl von Ossietzky Universitaet Oldenburg
* Ammerlaender Heerstr. 114-118
* 26111 Oldenburg, Germany
*
* CERIUUM
* Laboratoire SIM
* Faculte de Medecine
* 2 Avenue du Pr. Leon Bernard
* 35043 Rennes Cedex, France
*
* for CEN/TC251/WG4 as a contribution to the Radiological
* Society of North America (RSNA) 1993 Digital Imaging and
* Communications in Medicine (DICOM) Demonstration.
*
* THIS SOFTWARE IS MADE AVAILABLE, AS IS, AND NEITHER OFFIS,
* OLDENBURG UNIVERSITY NOR CERIUUM MAKE ANY WARRANTY REGARDING
* THE SOFTWARE, ITS PERFORMANCE, ITS MERCHANTABILITY OR
* FITNESS FOR ANY PARTICULAR USE, FREEDOM FROM ANY COMPUTER
* DISEASES OR ITS CONFORMITY TO ANY SPECIFICATION. THE
* ENTIRE RISK AS TO QUALITY AND PERFORMANCE OF THE SOFTWARE
* IS WITH THE USER.
*
* Copyright of the software and supporting documentation
* is, unless otherwise stated, jointly owned by OFFIS,
* Oldenburg University and CERIUUM and free access is hereby
* granted as a license to use this software, copy this
* software and prepare derivative works based upon this
* software. However, any distribution of this software
* source code or supporting documentation or derivative
* works (source code and supporting documentation) must
* include the three paragraphs of this copyright notice.
*
*/

```

Some other parts of this software within the dcmtk/dcmnet sub-package

related to the DICOM Upper Layer Protocol are derived from software developed for the RSNA'93 DICOM demonstration and kindly made available to us by the Mallinckrodt Institute of Radiology. Such software can be identified by the following copyright notice located in each affected source file:

```
/*
 * Copyright (C) 1993, RSNA and Washington University
 *
 * The software and supporting documentation for the Radiological
 * Society of North America (RSNA) 1993 Digital Imaging and
 * Communications in Medicine (DICOM) Demonstration were developed
 * at the
 *     Electronic Radiology Laboratory
 *     Mallinckrodt Institute of Radiology
 *     Washington University School of Medicine
 *     510 S. Kingshighway Blvd.
 *     St. Louis, MO 63110
 * as part of the 1993 DICOM Central Test Node project for, and
 * under contract with, the Radiological Society of North America.
 *
 * THIS SOFTWARE IS MADE AVAILABLE, AS IS, AND NEITHER RSNA NOR
 * WASHINGTON UNIVERSITY MAKE ANY WARRANTY ABOUT THE SOFTWARE, ITS
 * PERFORMANCE, ITS MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR
 * USE, FREEDOM FROM ANY COMPUTER DISEASES OR ITS CONFORMITY TO ANY
 * SPECIFICATION. THE ENTIRE RISK AS TO QUALITY AND PERFORMANCE OF
 * THE SOFTWARE IS WITH THE USER.
 *
 * Copyright of the software and supporting documentation is
 * jointly owned by RSNA and Washington University, and free access
 * is hereby granted as a license to use this software, copy this
 * software and prepare derivative works based upon this software.
 * However, any distribution of this software source code or
 * supporting documentation or derivative works (source code and
 * supporting documentation) must include the three paragraphs of
 * the copyright notice.
 */
```

The dcmjpeg sub-package includes an adapted version of the Independent JPEG Group Toolkit Version 6b, which is contained in dcmjpeg/libijg8, dcmjpeg/libijg12 and dcmjpeg/libijg16. This toolkit is covered by the following copyright. The original README file for the Independent JPEG Group Toolkit is located in dcmjpeg/docs/ijg_readme.txt.

```
/*
 * The authors make NO WARRANTY or representation, either express or implied,
 * with respect to this software, its quality, accuracy, merchantability, or
 * fitness for a particular purpose. This software is provided "AS IS", and you,
 * its user, assume the entire risk as to its quality and accuracy.
 *
 * This software is copyright (C) 1991-1998, Thomas G. Lane.
 * All Rights Reserved except as specified below.
 *
 * Permission is hereby granted to use, copy, modify, and distribute this
 * software (or portions thereof) for any purpose, without fee, subject to these
 * conditions:
 * (1) If any part of the source code for this software is distributed, then this
 * README file must be included, with this copyright and no-warranty notice
 * unaltered; and any additions, deletions, or changes to the original files
 * must be clearly indicated in accompanying documentation.
 * (2) If only executable code is distributed, then the accompanying
 * documentation must state that "this software is based in part on the work of
 * the Independent JPEG Group".
 * (3) Permission for use of this software is granted only if the user accepts
 * full responsibility for any undesirable consequences; the authors accept
 * NO LIABILITY for damages of any kind.
 *
 * These conditions apply to any software derived from or based on the IJG code,
 * not just to the unmodified library. If you use our work, you ought to
 * acknowledge us.
 *
 * Permission is NOT granted for the use of any IJG author's name or company name
 * in advertising or publicity relating to this software or products derived from
 * it. This software may be referred to only as "the Independent JPEG Group's
 * software".
 *
 * We specifically permit and encourage the use of this software as the basis of
```

```
* commercial products, provided that all warranty or liability claims are
* assumed by the product vendor.
*/
```

The code for the interpolatePixel() image scaling algorithm in module dcmimgle has been derived from code written by Jef Poskanzer for the "Extended Portable Bitmap Toolkit" (pbmplus10dec91) which has the following copyright:

```
/*
 * Copyright (C) 1991 by Jef Poskanzer.
 *
 * Permission to use, copy, modify, and distribute this software and its
 * documentation for any purpose and without fee is hereby granted, provided
 * that the above copyright notice appear in all copies and that both that
 * copyright notice and this permission notice appear in supporting
 * documentation. This software is provided "as is" without express or
 * implied warranty.
*/
```

The color quantization code in module dcmimage (dcmquant and the related classes) has been derived from code written by Jef Poskanzer for the NetPBM toolkit which has the following copyright:

```
/*
 * Copyright (C) 1989, 1991 by Jef Poskanzer.
 *
 * Permission to use, copy, modify, and distribute this software and its
 * documentation for any purpose and without fee is hereby granted, provided
 * that the above copyright notice appear in all copies and that both that
 * copyright notice and this permission notice appear in supporting
 * documentation. This software is provided "as is" without express or
 * implied warranty.
*/
```

The code for the OFStandard::strncpy and OFStandard::strlcat helper functions in ofstd/libsrc/ofstd.cc has been derived from the BSD implementation of strncpy() and strlcat() and which carries the following copyright notice:

```
/*
 * Copyright (c) 1998 Todd C. Miller <Todd.Miller(at)courtesan.com>
 * All rights reserved.
 *
 * Redistribution and use in source and binary forms, with or without
 * modification, are permitted provided that the following conditions
 * are met:
 * 1. Redistributions of source code must retain the above copyright
 * notice, this list of conditions and the following disclaimer.
 * 2. Redistributions in binary form must reproduce the above copyright
 * notice, this list of conditions and the following disclaimer in the
 * documentation and/or other materials provided with the distribution.
 * 3. The name of the author may not be used to endorse or promote products
 * derived from this software without specific prior written permission.
 *
 * THIS SOFTWARE IS PROVIDED ``AS IS'' AND ANY EXPRESS OR IMPLIED WARRANTIES,
 * INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY
 * AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL
 * THE AUTHOR BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL,
 * EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO,
 * PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS;
 * OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY,
 * WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR
 * OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF
 * ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.
*/
```

The code for the OFStandard::atof helper function in ofstd/libsrc/ofstd.cc has been derived from an implementation which carries the following copyright notice:

```
/*
 * Copyright 1988 Regents of the University of California
 * Permission to use, copy, modify, and distribute this software and
 * its documentation for any purpose and without fee is hereby granted,
 * provided that the above copyright notice appear in all copies. The
 * University of California makes no representations about the
```

```

* suitability of this software for any purpose. It is provided "as
* is" without express or implied warranty.
*
*
* The code for OFStandard::ftoa has been derived
* from an implementation which carries the following copyright notice:
*
* Copyright (c) 1988 Regents of the University of California.
* All rights reserved.
*
* Redistribution and use in source and binary forms are permitted
* provided that the above copyright notice and this paragraph are
* duplicated in all such forms and that any documentation,
* advertising materials, and other materials related to such
* distribution and use acknowledge that the software was developed
* by the University of California, Berkeley. The name of the
* University may not be used to endorse or promote products derived
* from this software without specific prior written permission.
* THIS SOFTWARE IS PROVIDED ``AS IS'' AND WITHOUT ANY EXPRESS OR
* IMPLIED WARRANTIES, INCLUDING, WITHOUT LIMITATION, THE IMPLIED
* WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.
*/

```

The "Base64" encoder/decoder in ofstd/libsrc/ofstd.cc has been derived from an implementation which carries the following copyright notice:

```

/*
* Copyright (c) 1999, Bob Withers - bwit(at)pobox.com
*
* This code may be freely used for any purpose, either personal or commercial,
* provided the authors copyright notice remains intact.
*/

```

The oflog sub-package is based on the log4cplus library which is covered by the following two copyright notices (for details see oflog/docs/LICENSE):

```

/*
* Copyright (C) 1999-2009 Contributors to log4cplus project.
* All rights reserved.
*
* Redistribution and use in source and binary forms, with or without modifica-
* tion, are permitted provided that the following conditions are met:
*
* 1. Redistributions of source code must retain the above copyright notice,
* this list of conditions and the following disclaimer.
*
* 2. Redistributions in binary form must reproduce the above copyright notice,
* this list of conditions and the following disclaimer in the documentation
* and/or other materials provided with the distribution.
*
* THIS SOFTWARE IS PROVIDED ``AS IS'' AND ANY EXPRESSED OR IMPLIED WARRANTIES,
* INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND
* FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE
* APACHE SOFTWARE FOUNDATION OR ITS CONTRIBUTORS BE LIABLE FOR ANY DIRECT,
* INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLU-
* DING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS
* OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON
* ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT
* (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF
* THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.
*/

```

```

// Module: Log4CPLUS
// File: appender.h
// Created: 6/2001
// Author: Tad E. Smith
//
//
// Copyright 2001-2009 Tad E. Smith
//
// Licensed under the Apache License, Version 2.0 (the "License");
// you may not use this file except in compliance with the License.
// You may obtain a copy of the License at
//
// http://www.apache.org/licenses/LICENSE-2.0
//

```

```
// Unless required by applicable law or agreed to in writing, software
// distributed under the License is distributed on an "AS IS" BASIS,
// WITHOUT WARRANTIES OR CONDITIONS OF ANY KIND, either express or implied.
// See the License for the specific language governing permissions and
// limitations under the License.
```

The dcmjpls sub-package is based on the CharLS library, which is contained in dcmjpls/libcharls. This toolkit is covered by the following copyright:

```
/*
 * The CharLS library is available under the following license:
 * -----
 * Copyright (c) 2007-2010, Jan de Vaan
 * All rights reserved.
 *
 * Redistribution and use in source and binary forms, with or without
 * modification, are permitted provided that the following conditions are met:
 *
 * * Redistributions of source code must retain the above copyright notice, this
 *   list of conditions and the following disclaimer.
 *
 * * Redistributions in binary form must reproduce the above copyright notice,
 *   this list of conditions and the following disclaimer in the documentation
 *   and/or other materials provided with the distribution.
 *
 * * Neither the name of my employer, nor the names of its contributors may be
 *   used to endorse or promote products derived from this software without
 *   specific prior written permission.
 *
 * THIS SOFTWARE IS PROVIDED BY THE COPYRIGHT HOLDERS AND CONTRIBUTORS "AS IS"
 * AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE
 * IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE
 * DISCLAIMED. IN NO EVENT SHALL THE COPYRIGHT OWNER OR CONTRIBUTORS BE LIABLE FOR
 * ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES
 * (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES;
 * LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON
 * ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT
 * (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS
 * SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.
 */
```